


How to change tube for your hearing aid (Mould)

Tubing will need to be replaced every 4-6 months and will become discoloured, hard/ brittle over time.

Separate the mould from your hearing aid

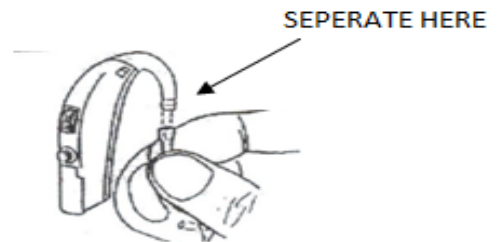

Remove the old tubing from the ear mould by pulling out with pliers. Save the old piece of tubing until step 7.




3. Cut a long sharp point at one end of the tubing with scissors.




Thread the pointed end of the tubing through the hole in the ear mould (if there are 2 holes in the ear mould place the tubing through the hole at the top)



Pull the tubing through with pliers until the non-cut end of the tubing points upwards and over the top of the ear mould. (See diagram)

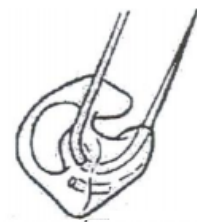



Cut the pointed end off flush where it meets with the end of the ear mould (take care not to cut the ear mould)

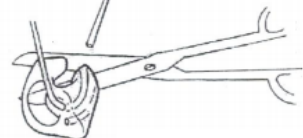

7. Use the old piece of tubing saved from step 2 to measure the correct length (see diagram)



8. reattach the ear mould to your hearing aid.

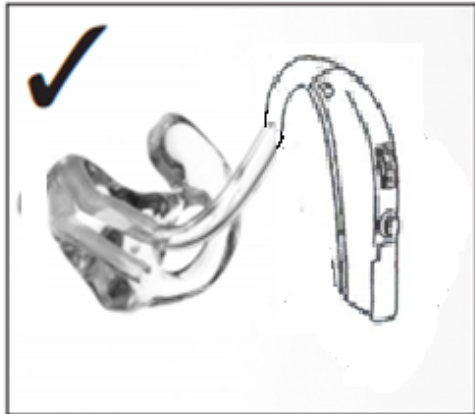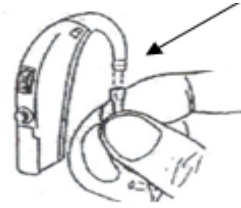

9. The mould should face away from the hearing aid and not towards as shown in the above diagrams.

If you have any problems please contact us on the details below or alternatively drop into our walk- in service either on a Wednesday or a Friday between 9.00am - 12.00 noon.

**Audiology
Modality Community Services**

Address: Attwood Green Medical Centre, 15 Bath Row, Birmingham, B15
1LZ

Website: www.modalitycommunityservices.com

Phone 0121 250 1592 opt 5

Email: modality.audiology@nhs.net